

COMPACT T Short Notch Automatic Soft Tissue Biopsy Device

One-handed design simplifies loading and operation whilst allowing for convenient manipulation of the ultrasound transducer.

- Programmable: releases the stylet allowing time for tissue to settle in notch (ideal for denser tissue or when there is a desire to verify notch position) and activates the cutting cannula
- Sample notch 10mm for critical procedures
- Thin-wall cannula for larger core
- Colour coded for quick gauge size identification
- Depth markings etched on cannula for precise positioning
- Loaded status indicator
- Plastic safety protection to avoid accidental firing




COMPACT

PRODUCT CODE:	SIZE
VCT11410	14G x 10cm
VCT11413	14G x 13cm
VCT11415	14G x 15cm
VCT11810	18G x 10cm
VCT12010	20G x 10cm

(Box of x10 pieces)

Needles with different diameters and lengths may be provided on request.

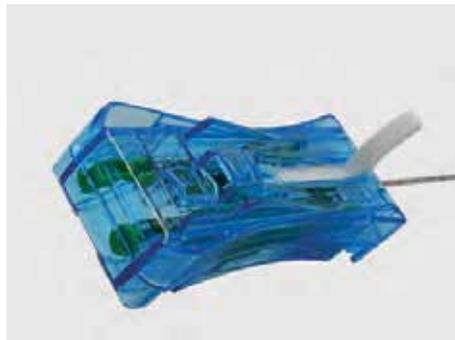

Colour code and plastic safety protection