


BIOMARC IGRT Fiducial Marker for Radiation Therapy

- Used to target localisation for lesions or organs which move in respect to external anatomy
- Pure gold cylinders version specially knurled to inhibit migration
- Pyrolytic carbon coated zirconium oxide or tantalum barbell version to reduce the artifacts
- Provide visualisation under fluoroscopy, X-Ray, CT, MRI, ultrasound and other imaging modalities
- Preloaded in a pre-waxed tipped needle or in a wax-free special version
- Reduces the artifact in the MV setting (BiomarC EX)
- Ideal for proton or photon therapy (BiomarC Pro)



BioMarC GOLD

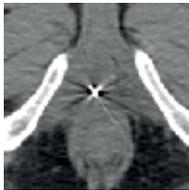

BioMarC KV ENH



BioMarC KV




BioMarC EX



BioMarC PRO



BioMarC SECURE

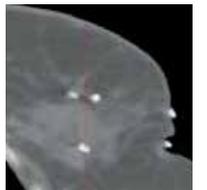

BIOMARC IGRT

PRODUCT CODE:	SIZE	DESCRIPTION	QTY
040116	17G x 20cm	1.2 x 3mm Gold Marker	3
040118	17G x 20cm	1.2 x 4.2mm Gold Marker – Ergonomic handle	10
040141	18G x 20cm	0.9 x 3mm Gold Marker	3
040125	17G x 30cm	1.2 x 3mm Gold Marker	3
040117	18G x 20cm	1 x 5mm Carbon KV ENH Marker	3
040113	18G x 20cm	1 x 3mm Carbon KV Marker	3
040146	18G x 20cm	1 x 5mm Carbon EX Marker	3
040148	18G x 20cm	1 x 10mm Carbon PRO Marker	1
040150	Sutured	1.7 x 4.6mm Carbon Secure Marker	10

Needles with different diameters and lengths may be provided on request.